OTICON | Zircon

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miniRITE R

60 65 100 105

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
miniRITE R

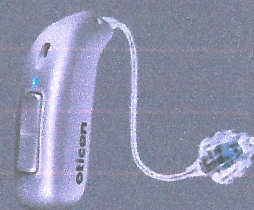60 85 100 105

Oticon Zircon miniRITE R offers a discreet design powered by a rechargeable lithium-ion battery. The style features telecoil, and a double push-button. It is a Made for iPhone® hearing aid and compatible with the new Android protocol for Audio Streaming for Hearing Aids (ASHA) - making it possible to stream directly from iPhone, iPad®, iPod touch® and selected Android™ devices.

OpenSound Navigator™ provides access to speech in 360° making the listener more easily aware of what is going on in the surroundings.

Speech Guard™ provides more natural and clear speech sounds making the details in speech stand out more.

The Polaris™ platform provides a tremendous speed and memory capacity for audiological processing and connectivity options. New features can be added and updates performed wirelessly.

General features:

- Digital Programmable
- Automatic or Manual Volume Control
- Maximum OutputControl System
- MPO-Maximum Power Output
- GC-Gain Control
- AGC-Automatic Gain Control
- Noise Reduction
- Feedback Management
- Dual Microphone
- FM Compatible (with Telecoil)
- 4 Programs

	Zircon 1	Zircon 2
Speech Understanding	OpenSound Navigator™	•
	- Balancing power effect	40%
	- Max. noise removal difficult/simple	6 dB / 0 dB
	Multiband Adaptive Directionality	-
	Noise Reduction	-
	Speech Guard™	•
	Single Compression	-
	Frequency lowering	Speech Rescue™
Sound Quality	Fitting Bandwidth*	8 kHz
	Bass Boost (streaming)	•
	Processing Channels	48
Listening Comfort	Feedback Management	SuperShield & Feedback shield
	Transient Noise Management	On/Off
	Wind Noise Management	•
Personalisation & Optimising Fitting	Fitting Bands	14
	Multiple Directionality options	•
	Adaptation Management	•
	Oticon Firmware Updater	•
	Fitting Formulas	NAL-NL1/NAL-NL2, DSL 5.0
Connecting to the world	Hands-free communication**	•
	Direct streaming***	•
	Oticon ON app & Oticon RemoteCare app	•
	ConnectClip	•
	EduMic	•
	Remote Control 3.0	•
	TV Adapter 3.0	•
	Phone Adapter 2.0	•
	Tinnitus SoundSupport™	•
	CROS/BiCROS support	•

*Bandwidth accessible for gain adjustments during fitting

**Available for Oticon Zircon from FW 1.1 with selected iPhone models

***From iPhone®, iPad®, iPod touch®, and selected Android™ devices

Operating and charging conditions

Temperature: +5°C to +40°C (41°F to 104°F)
Relative humidity: 5% to 93%, non-condensing
Atmospheric pressure: 700 hPa to 1060 hPa

Storage and transportation conditions

Temperature and humidity should not exceed the below limits for extended periods during transportation and storage.

Transport

Temperature: -20°C to +60°C (-4°F to 140°F)
Relative humidity: 5% to 93%, non-condensing
Atmospheric pressure: 700 hPa to 1060 hPa

Storage

Temperature: -20°C to +30°C (-4°F to 86°F)
Relative humidity: 5% to 93%, non-condensing
Atmospheric pressure: 700 hPa to 1060 hPa

Apple, the Apple logo, iPhone, iPad, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries.

Made for
iPhone | iPad | iPod

Works with
android

IP68

For information on compatibility, please visit www.oticon.global/compatibility

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

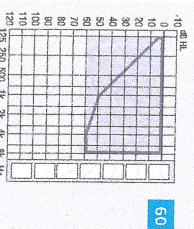
ZCC Coupler

Measured according to
ANSI S3.22 2014, IEC 60118-0:2015
and IEC 60318-5:2008

☐ Mould, BASS & Power dome
☐ Operates dome

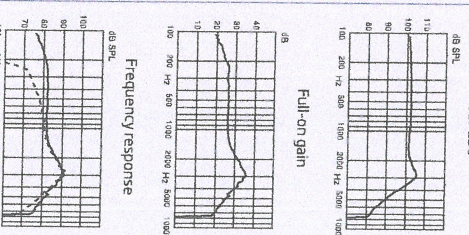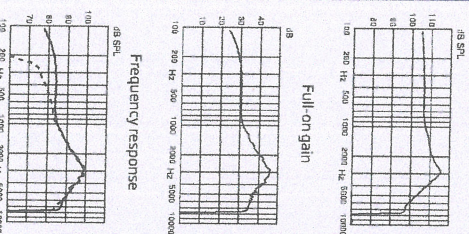Technical information
Omnidirectional mode is used unless otherwise stated.

Acoustic input: 60 dB SPL
Magnetic input: 31.5 mA/m

Acoustic input: 60 dB SPL
Magnetic input: 31.5 mA/m

OSPL90

Peak
1600 Hz
HFA-OSPL90

116 dB SPL
110 dB SPL
110 dB SPL

106 dB SPL
102 dB SPL
103 dB SPL

Full-on gain¹⁾

Peak
1600 Hz
HFA-FOG

46 dB
37 dB
38 dB

36 dB
29 dB
30 dB

Reference test gain

100-7500 Hz

31 dB

26 dB

Frequency range

100-7500 Hz

68 dB SPL
88 dB SPL

83/83 dB SPL

Telecoil output (1600 Hz)

1 mA/m field
10 mA/m field

500 Hz
800 Hz
1600 Hz

83/83 dB SPL

Total harmonic distortion (Input 70 dB SPL)

500 Hz
800 Hz
1600 Hz

< 2 %
< 3 %
< 2 %

< 2 %
< 2 %
< 2 %

Equivalent input noise level

Omni
Dir

19 dB SPL
26 dB SPL

17 dB SPL
29 dB SPL

Battery

Lithium-ion

24

Expected operating time, hours²⁾

Lithium-ion

24

KOBI BETA
TRAKIMAN İNHALECE DEN TÜRKİYEYE GELİMLİSTİR
EMRAH BÜKE YENİMLİ MÜTEREM

Oticon Zircon 2

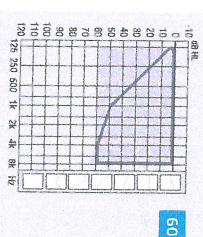
miniRITE R60

Ear Simulator

Measured according to
IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015,
IEC 60118-1:1985+AMD1:1998/CS and
IEC 60318-4:2010



ZCC Coupler

Measured according to
ANSI S3.22 2014, IEC 60118-0:2015
and IEC 60318-5:2008

☐ Mould, BASS & Power dome
☐ Operates dome

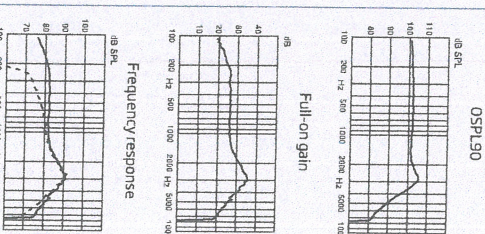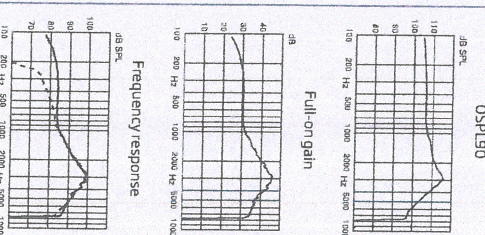Technical information
Omnidirectional mode is used unless otherwise stated.

Acoustic input: 60 dB SPL
Magnetic input: 31.5 mA/m

Acoustic input: 60 dB SPL
Magnetic input: 31.5 mA/m

OSPL90

Peak
1600 Hz
HFA-OSPL90

116 dB SPL
110 dB SPL
110 dB SPL

106 dB SPL
102 dB SPL
103 dB SPL

Full-on gain¹⁾

Peak
1600 Hz
HFA-FOG

46 dB
37 dB
38 dB

36 dB
29 dB
30 dB

Reference test gain

100-7500 Hz

31 dB

26 dB

Frequency range

100-7500 Hz

68 dB SPL
88 dB SPL

83/83 dB SPL

Telecoil output (1600 Hz)

1 mA/m field
10 mA/m field

500 Hz
800 Hz
1600 Hz

83/83 dB SPL

Total harmonic distortion (Input 70 dB SPL)

500 Hz
800 Hz
1600 Hz

< 2 %
< 3 %
< 2 %

< 2 %
< 2 %
< 2 %

Equivalent input noise level

Omni
Dir

19 dB SPL
26 dB SPL

17 dB SPL
29 dB SPL

Battery

Lithium-ion

24

Expected operating time, hours²⁾

Lithium-ion

24

¹⁾ Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+AMD1:1994 but without influence of feedback.

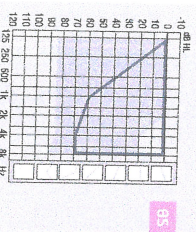
²⁾ Expected operating time for rechargeable battery depends on use pattern, active feature set, hearing loss, sound environment, battery age and use of wireless accessories.

¹⁾ Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+AMD1:1994 but without influence of feedback.

²⁾ Expected operating time for rechargeable battery depends on use pattern, active feature set, hearing loss, sound environment, battery age and use of wireless accessories.

Oticon Zircon 1

miniRITE R 85


☐ Mold, Batts & Power dome
☐ Operates dome


Technical information
Omnidirectional mode is used unless otherwise stated.

Ear Simulator
Measured according to
IEC 60118-0:1983/AM1:1994, IEC 60119-0:2015,
IEC 60118-1:1995-AM01:1998 CSV and
IEC 60318-4:2010

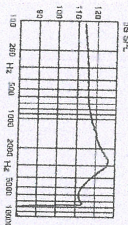OSPL90


Full-on gain

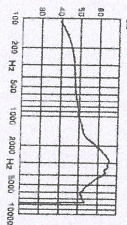
Frequency response

ZCC Coupler
Measured according to
ANSI S3.22:2016, IEC 60118-0:2015
and IEC 60318-5:2008

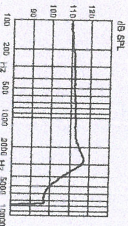OSPL90

Full-on gain

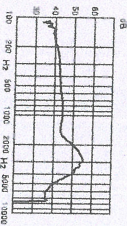Frequency response

OSPL90

Peak
1600 Hz
HFA-OSPL90
122 dB SPL

Full-on gain¹
Peak
1600 Hz
HFA-FOG
56 dB

Reference test gain
Frequency range

Telecoil output (1600 Hz)

1 mA/m field
10 mA/m field
SPLITS L/R
500 Hz
800 Hz
1600 Hz
Omni
Dir

Total harmonic distortion (Input 70 dB SPL)

Equivalent input noise level

Battery
Expected operating time, hours²

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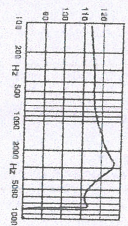
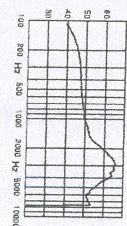
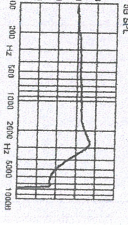
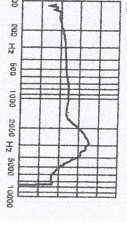
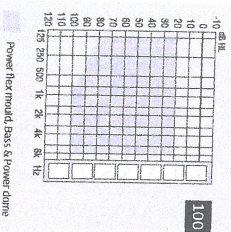
1) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+AM1:1994 but without influence of feedback.
2) Expected operating time for rechargeable battery depends on use pattern, active feature set, hearing loss, sound environment, battery age and use of wireless accessories.

Oticon Zircon 1

miniRITE R100



Ear Simulator
Measured according to
IEC 60118-0:1983/AM01:1994 and IEC 60118-0:2015,
IEC 60118-1:1995+AM01:1999 CSV and
IEC 60318-4:2010

ZCC Coupler
Measured according to
ANSI S3.22:2014, IEC 60118-0:2015
and IEC 60318-5:2008

OSPL90

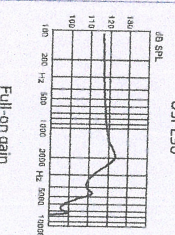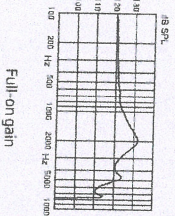OSPL90



Technical information
Omnidirectional mode is used unless otherwise stated.
Warning to the hearing aid dispenser
The maximum output capability of the hearing aid may exceed 132 dB SPL (IEC 711). Special care should be exercised in selecting and fitting the hearing aid, as there may be risk of impairing the remaining hearing of the hearing aid user.

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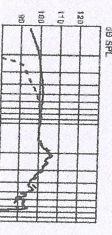
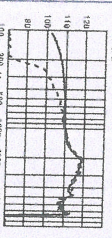
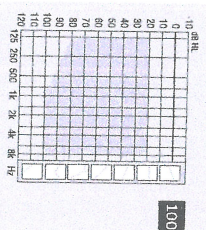
ZCC Coupler
Measured according to
ANSI S3.22:2014, IEC 60118-0:2015
and IEC 60318-5:2008

OSPL90

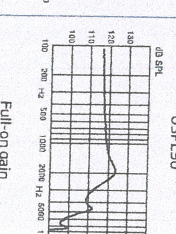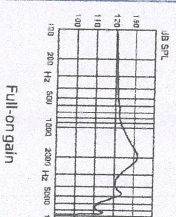OSPL90



Technical information
Omnidirectional mode is used unless otherwise stated.
Warning to the hearing aid dispenser
The maximum output capability of the hearing aid may exceed 132 dB SPL (IEC 711). Special care should be exercised in selecting and fitting the hearing aid, as there may be risk of impairing the remaining hearing of the hearing aid user.

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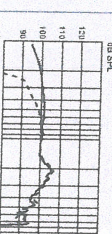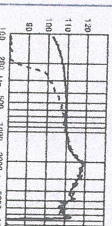Peak
1600 Hz
HFA-FOG

24

1) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+AM1:1994 but without influence of feedback.

2) Expected operating time for rechargeable battery depends on use pattern, active features set, hearing loss, sound environment, battery age and use of wireless accessories.

1) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+AM1:1994 but without influence of feedback.

2) Expected operating time for rechargeable battery depends on use pattern, active features set, hearing loss, sound environment, battery age and use of wireless accessories.

Oticon Zircon 1

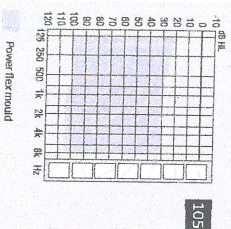
miniRITE R105

Ear Simulator

Measured according to
IEC 60318-01:1985/AMD1:1994, IEC 60318-0:2015,
IEC 60318-1:1995/AMD1:1998 CS and
IEC 60318-4:2010

ZCC Coupler

Measured according to
ANSI S3.22:2014, IEC 60318-0:2015
and IEC 60318-5:2006

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

Warning to the hearing aid dispenser

The maximum output capability of the hearing aid may exceed 132 dB SPL (IEC 711). Special care should be exercised in selecting and fitting the hearing aid, as there may be risk of impairing the remaining hearing of the hearing aid user.

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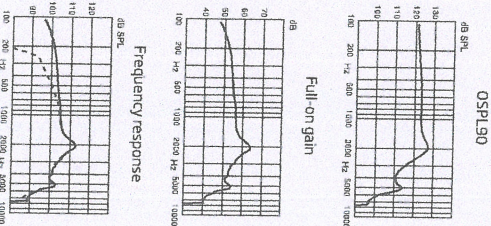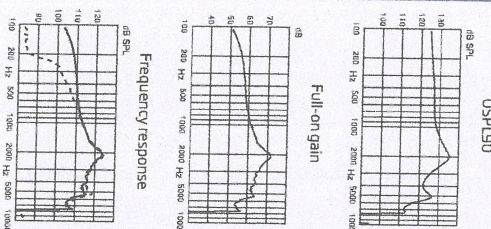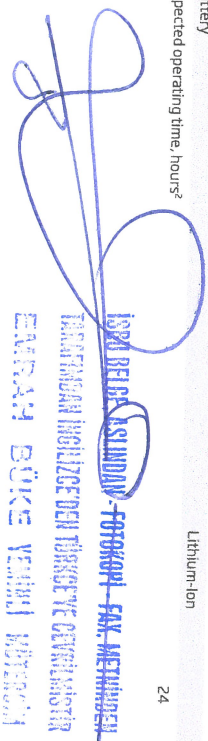
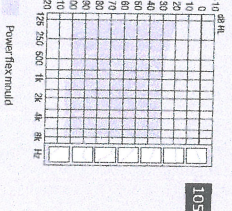
ZCC Coupler

Measured according to
ANSI S3.22:2014, IEC 60318-0:2015
and IEC 60318-5:2006

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

Warning to the hearing aid dispenser

The maximum output capability of the hearing aid may exceed 132 dB SPL (IEC 711). Special care should be exercised in selecting and fitting the hearing aid, as there may be risk of impairing the remaining hearing of the hearing aid user.

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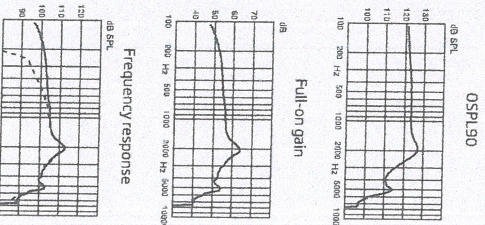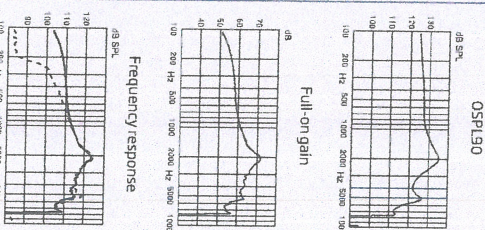Expected operating time, hours²

Peak
1600 Hz
HFA-FOG

¹) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60318-01:1985/AMD1:1994 but without influence of feedback.

²) Expected operating time for rechargeable battery depends on use pattern, active features set, hearing loss, sound environment, battery age and use of wireless accessories.

¹) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60318-01:1985/AMD1:1994 but without influence of feedback.

²) Expected operating time for rechargeable battery depends on use pattern, active features set, hearing loss, sound environment, battery age and use of wireless accessories.

İSRAİL BELEDİYESİNDEN FOTOKOPİ - FAX MEYDANCI
TARAFINDAN İNCELENİP TIRAGEYA ÇEVİRİLMİŞTİR
EMİNE ÇİĞDEM YILMAZ / 0532 304 18 20 00 04

Headquarters
Oticon A/S
Kongebakken 9
DK-2765 Smørum
Denmark

SBO Hearing A/S
Kongebakken 9
DK-2765 Smørum
Denmark

OTICON A/S
KONGEBAKKEN 9
DK-2765 SMØRUM
DENMARK

247593TR-UK / 2022.01.03 / v1

www.oticon.global

Oticon is part of the Demant Group.

oticon
life-changing technology